

# Rotary vane pump RZ 6 +FO +VS 16 +Set DCP+VSP 3000

[Link to the website](#)

## Product description

These powerful rotary vane pumps feature a particularly compact design and low weight for pumps of this capacity. They are the ideal solution for a wide range of laboratory and process applications that require low ultimate vacuum at medium to increased gas flow rate. Various packages including pump, oil mist filter, etc. are available.

## Performance features

- high flow rates even at vacuum levels approaching ultimate vacuum
- high water vapor tolerance due to efficient gas ballast; very good ultimate vacuum even with gas ballast
- vacuum-tight at switch-off; external anti-suckback valve not needed
- large oil volume: Long intervals between oil changes
- ease of maintenance due to telescopic design

Ordering data

| Product name                      | Nominal mains voltage   mains frequency | Plug type | Articleno. |
|-----------------------------------|-----------------------------------------|-----------|------------|
| RZ 6 +FO +VS 16 +Set DCP+VSP 3000 | 230 V   50-60 Hz                        | CEE       | 20798150   |
| RZ 6 +FO +VS 16 +Set DCP+VSP 3000 | 230 V   50-60 Hz                        | CH, CN    | 20798151   |

Items supplied

Pump oil filled and completely mounted, ready for use, with manual, with oil mist filter FO, valve VS 16 and Pirani vacuum gauge DCP 3000 +VSP3000


Accessories

- Rubber vacuum tubing DN 10 mm (20686002)
- Metal vacuum tubing with small flanges KF DN 16 (20673336)
- Separator inlet side AK for RE/RZ 5/6/2.5 (20698006)
- Maintenance kit RE 2.5 - RZ 6 (20648977)


## Technical data

| Max. pumping speed 50 Hz (mbar)                                                                                  | 5.7 m <sup>3</sup> /h                                    |
|------------------------------------------------------------------------------------------------------------------|----------------------------------------------------------|
| Max. pumping speed 60 Hz (cfm)                                                                                   | 4.0 cfm                                                  |
| Ultimate vacuum                                                                                                  | 2 x 10 <sup>-3</sup> mbar / 1.5 x 10 <sup>-3</sup> torr  |
| Ultimate vacuum with gas ballast                                                                                 | 1 x 10 <sup>-2</sup> mbar / 0.75 x 10 <sup>-2</sup> torr |
| Ultimate partial vacuum without gas ballast                                                                      | 4 x 10 <sup>-4</sup> mbar / 3 x 10 <sup>-4</sup> torr    |
| Water vapor tolerance                                                                                            | 40 mbar / 30 torr                                        |
| Number of stages                                                                                                 | 2                                                        |
| Vacuum controller                                                                                                | DCP 3000 / VSP 3000                                      |
| Permissible ambient temperature range (operation)                                                                | 12 - 40 °C                                               |
| Permissible ambient temperature range (storage)                                                                  | -10 - 60 °C                                              |
| Max. outlet pressure (abs.)                                                                                      | 1.1 bar                                                  |
| Oil capacity min./max.                                                                                           | 0.34 / 0.73 l                                            |
| Inlet connection                                                                                                 | Small flange KF DN 16 / hose nozzle DN 10 mm             |
| Outlet connection                                                                                                | Hose nozzle DN 8-10 mm                                   |
| Rated motor power                                                                                                | 0.18 kW                                                  |
| Rated speed 50/60 Hz                                                                                             | 1500/1800 min <sup>-1</sup>                              |
| Protection class IEC 60529                                                                                       | IP 40                                                    |
| Dimensions (L x W x H)                                                                                           | 420 mm x 230 mm x 240 mm                                 |
| Weight                                                                                                           | 17.8 kg                                                  |
| Noise (sound pressure level) in dBA at 50 Hz or 1500 rpm/62% (VARIO)/1500 rpm (VARIO-SP)/12500 Upm (VACUU·PURE®) | 50 dBA                                                   |
| NRTL certification                                                                                               | No                                                       |


## Pumping speed graph at 50 Hz

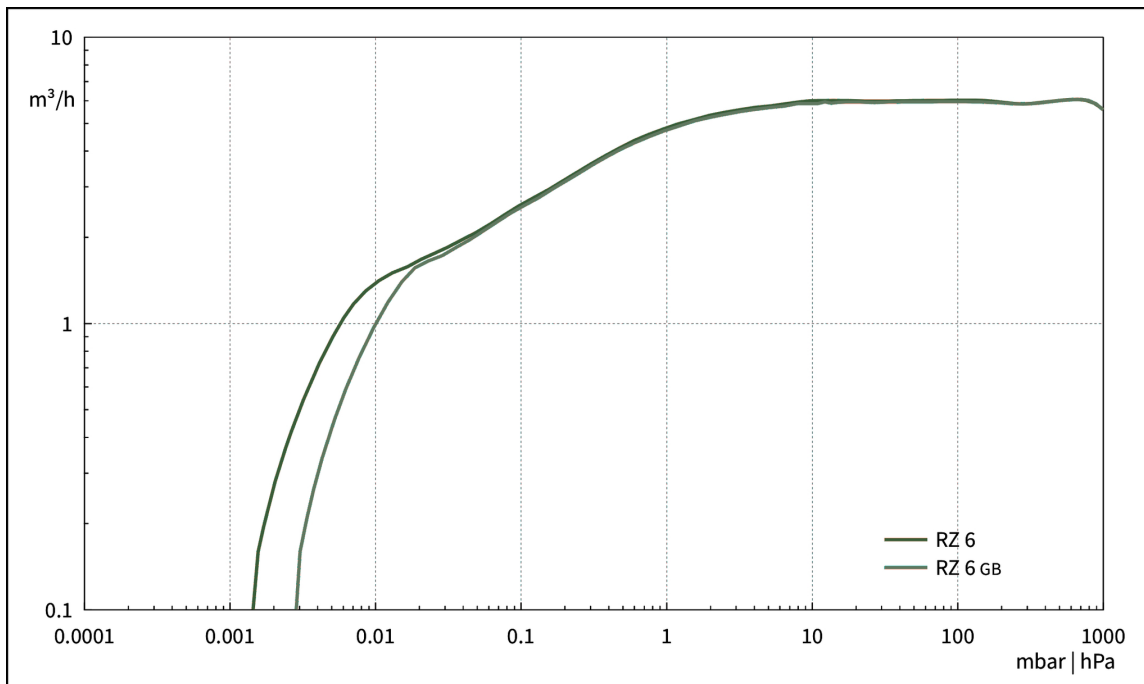


## Pump down graph at 50 Hz (10 l volume)

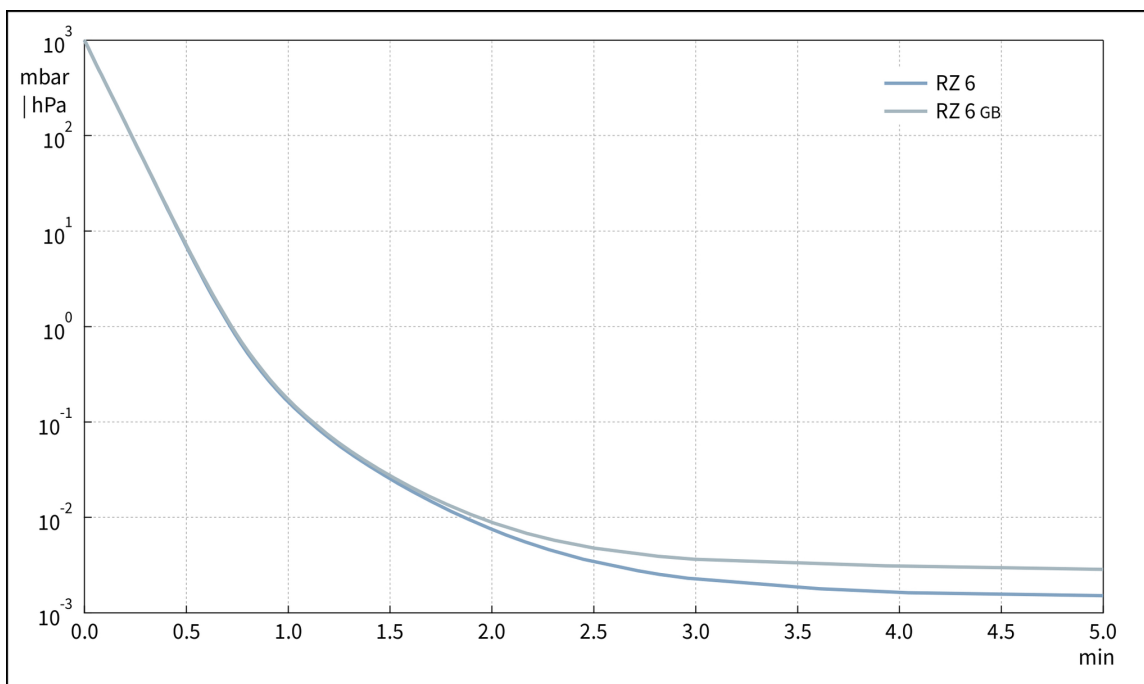


## Pumping speed graph at 60 Hz

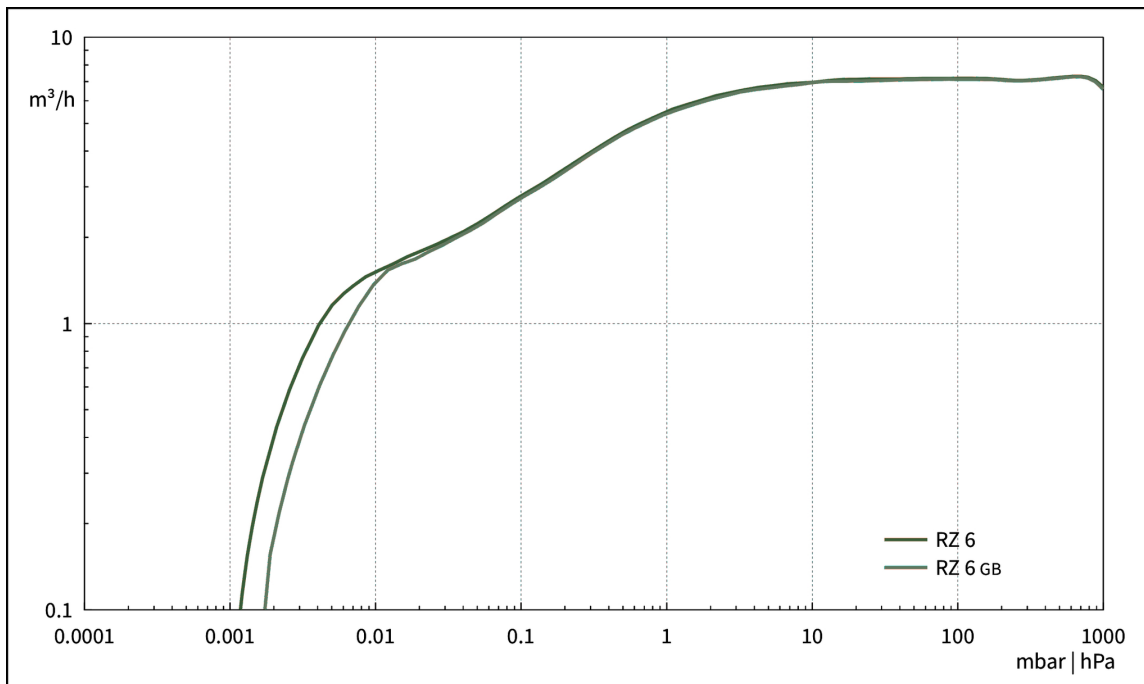

## Pump down graph at 60 Hz (10 l volume)

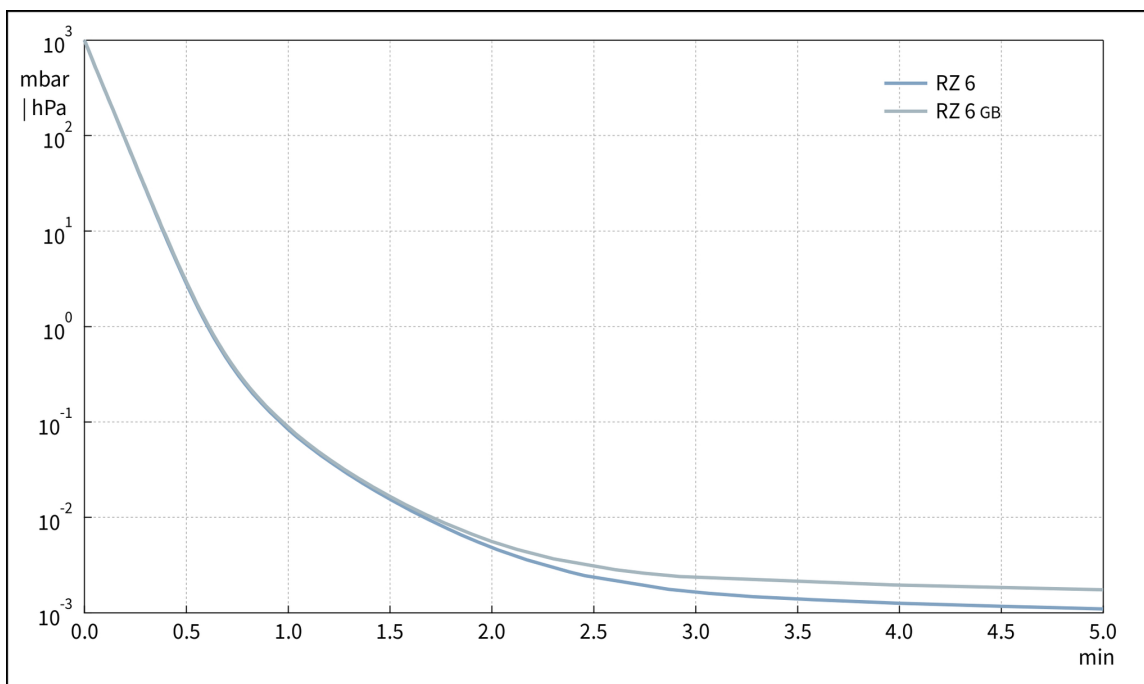

## Dimension sheet

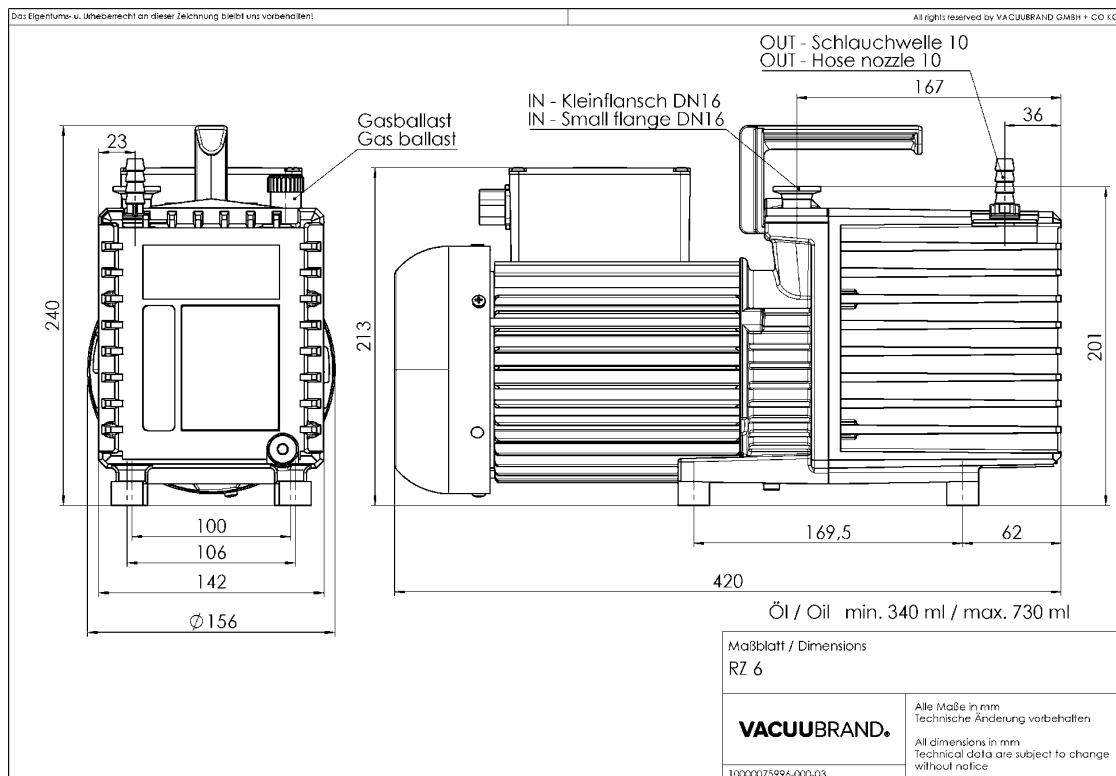

VACUUBRAND GMBH + CO KG  
P.O. Box 1664  
97877 Wertheim  
T +49 9342 808 5550  
F +49 9342 808 5555

info@vacuubrand.com  
www.vacuubrand.com

**VACUUBRAND®**

VACUUBRAND is part of the Brand group, a manufacturer of high-quality and innovative laboratory equipment and vacuum technology with a global team of approx. 1,000 employees.

With highly motivated and qualified employees, we manufacture the world's most comprehensive product family for vacuum generation, measurement and control for the rough and fine vacuum range in the laboratory.

**All with one goal:**

Enabling our customers to run reliable, safe and efficient processes. This is why we commit ourselves to continuous innovation, first-class performance, and outstanding quality.

Engineering excellence is central to our ability to deliver products and solutions renowned for their quality, durability, and reliability. Innovative technology, high quality materials and extensive testing deliver unmatched performance - engineered and manufactured in Germany.

Together with you, our experts will find the best solution for your application and provide you with first-class service along the way.

© VACUUBRAND GMBH + CO KG

BRAND (Shanghai) Trading Co., Ltd.  
Shanghai, China

T +86 21 6422 2318  
info@brand.com.cn  
china.brand.com.cn

BRAND Scientific Equipment Pvt. Ltd.  
Mumbai, India

T +91 22 42957790  
customersupport@brand.co.in  
www.brand.co.in

BRANDTECH® Scientific, Inc.  
Essex, CT. United States of America

T +1 860 767 2562  
info@brandtech.com  
www.brandtech.com